


# DISPOSABLE LIGATING CLIP

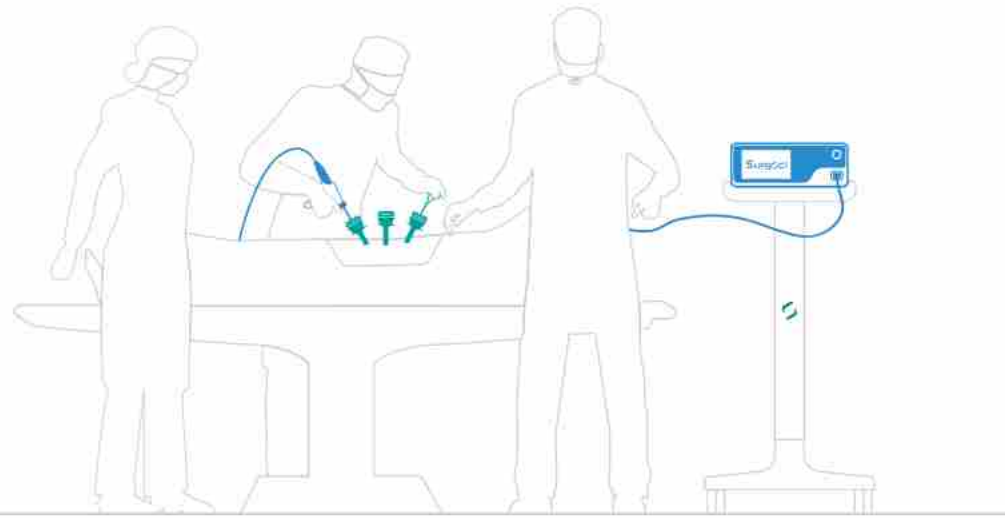

**Surgsci** SURGSCI MEDICAL LTD.

## SURGSCI MEDICAL LTD.

Tel: 0755-23707593 Email: [service@surgsci.com](mailto:service@surgsci.com)

### R&D Center:

Area A, 4/F, Summer Valley, Mei Sheng Hui Gu Technology Park, No. 83 Dabao Road, Shanghe Community, Xin'an Street, Bao'an District, Shenzhen, China

### Manufacturing Center:

301, No.3 Workshop, Run Heng Ding Feng Industrial Factory, #1 of Liuxian 3rd Road, Xin'an Street, Bao'an District, Shenzhen, Guangdong, China



contact us



whatsapp



*MEDICAL  
GRADE POM*

*High Level Biocompatibility and Safety*




DISPOSABLE LIGATING CLIP



**3** SIZES AVAILABLE FOR  
VARIOUS SURGERIES

**SURGSCI MEDICAL LTD.**

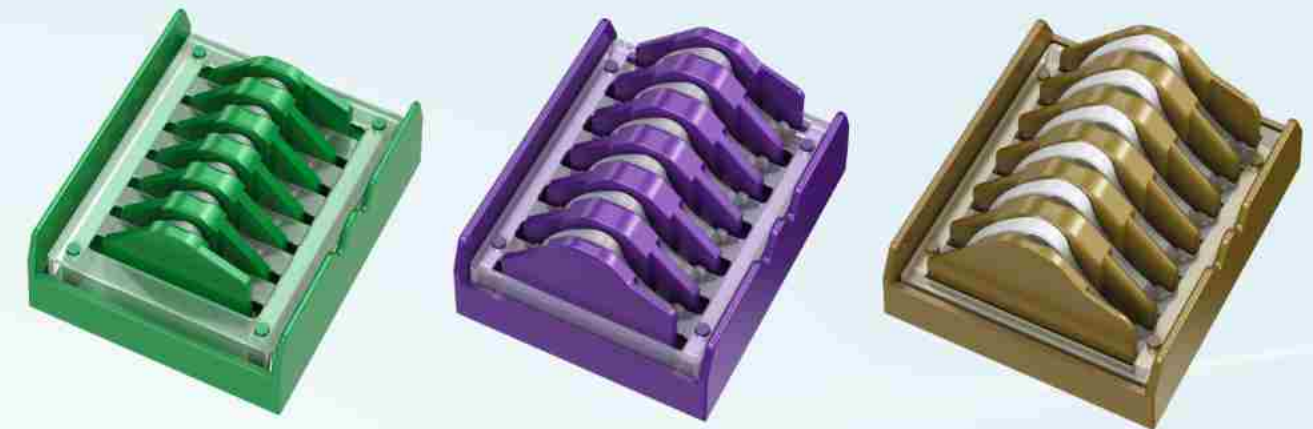

## SPECIFICATION OF DISPOSABLE LIGATING CLIP

	Medium-Large	Large	Extra-Large
			
<b>ML</b>	<b>L</b>	<b>XL</b>	
	Closed Length (formed): 3~5mm	Closed Length (formed): 4~8mm	Closed Length (formed): 9~12mm
	Suitable for laparoscopic cholecystectomy	General model, the most popular one on the market	Suitable for obstetrics, urology, and gastrointestinal surgery


## FEATURES AND ADVANTAGES

Disposable ligating clips are used to occlude blood vessels or other tubular structures during surgery and are used to replace traditional suture ligation. Compared with traditional surgical hemostats, the disposable ligating clip is smaller in size and more flexible in operation. It is more suitable for places that are difficult to reach with hemostats, such as the depth of narrow tissues.



### »»» Excellent Ligating Characteristics

Safety lock design with wide closed length. During the process of tissue detumescence, tension is always maintained, and it is not easy to loosen or slip. It is suitable for various operations

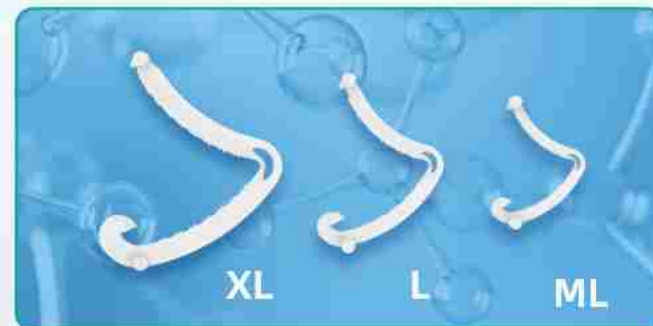







### »»» Medical POM Material

Made of medical-grade hydrolysis-resistant polyoxymethylene material, which has good biocompatibility and safety. It is radiolucent without imaging interference. This product is sterilized by ethylene oxide (EO) and is for one-time use.

## CLIP APPLIER & CLIP REMOVER

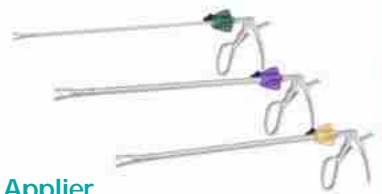Supplied with pistol grip clip applicator and clip remover. There are three sizes of clip applicators and two types of clip remover (as shown in the list below).

Clip applicators and removers are not sterile. The product needs to be sterilized before use and can be reused.

Product Type	Product Code	Shaft Diameter	Shaft Length
 Clip Applicator	 CA-03ML	5mm	330mm
	 CA-03L	10mm	330mm
	 CA-03XL	10mm	330mm
 Clip Remover	 CR-03ML	5mm	330mm
	 CR-03XL	10mm	330mm